Spin Air Basic

Instructions for Use

Air sampler

For units with SN 10005532/1065 and above



EN

Doc. Doc. 50006171-05

Revision date: October 2026

All Rights Reserved © 2026

Content

Indication of the Introduced Modifications **3**

1. Introduction **4**

1.1. Intended Use of the Spin Air Basic 4

1.2. Device Description 4

1.3. General Information 4

1.4. Requirements for Spin Air Basic Users 5

2. Safety Information **6**

2.1. Proper Use 6

2.2. Electrical Safety 7

2.3. Environment 8

2.4. Biological Safety 8

2.5. Chemicals 9

2.6. Maintenance Safety 9

2.7. Waste Disposal 10

2.8. Symbols on the Spin Air Basic Device 11

3. General Description **12**

3.1. Device Overview 12

3.2. Control Panel 13

3.3. Accessories of the Device 13

4. Installation **15**

4.1. Unpack the Device 15

4.2. Site Requirements 15

4.3. Power Cable Connection 15

5. Operating Procedures **16**

5.1. Sampling 16

5.2. Battery charging 17

5.3. Advanced Menu 18

6. Calibration **19**

7. Troubleshooting **20**

8. Maintenance **21**

8.1. Cleaning Procedure 21

APPENDIX A Ordering Codes **22**

APPENDIX B Technical Data **23**

APPENDIX C Warranty **25**

APPENDIX D Waste Electrical and Electronic Equipment (WEEE) **26**

APPENDIX E RoHS Statement **27**

APPENDIX F Statistical Corrections when Sampling Without Disk Rotation 28
APPENDIX G Table of Statistical Corrections for Sampling Using a 55/60 mm Dish 29
APPENDIX H Table of Statistical Corrections for Sampling Using a 90 mm Dish 31
APPENDIX I EU Declaration of Conformity 33
APPENDIX J Certificate of Compliance EN 17141:2020 34

Indication of the Introduced Modifications

If it's the first time that you use Spin Air Basic, please read the entire document carefully.

If it's not the first time that you use Spin Air Basic and you use it before, in the table below you can read the indication of the introduced modifications.

Revision	Description of change
05	Update format
	Update photo of front page
	Add SN number on front page
	Update sections 1.3.1, 3.1, 5.1, 5.2, 3.2

1. INTRODUCTION

Thank you for choosing the Spin Air Basic. We are confident that this instrument will become an integral part of your laboratory.

Before using the Spin Air Basic, it is essential that you read these Instructions for Use carefully. Following the instructions and safety information in this Instructions for Use will ensure safe operation and maintain the system in a safe condition.

1.1. Intended Use of the Spin Air Basic

The Spin Air Basic is intended to be used in the pharmaceutical, medical device, food, and cosmetic manufacturers, as well as healthcare facilities to assess airborne bacterial and fungal bioburden.

1.2. Device Description

The Spin Air Basic is a portable microbial air sampler that uses the unrivaled Spin technology. This sampling technology uses 100% of the Petri plate agar surface to plate microorganisms improving data statistical significance greatly and avoiding the use of colony count correction tables.

These samplers are compact and portable, they can be attached to a tripod for adequate sampling direction. A handy countdown function avoids operator interference with results.

The Spin Air Basic allows complying with USP 797 and 1116 regulations and is EN 17141:2020 Annex E certified, see APPENDIX J Certificate of Compliance EN 17141:2020.

1.3. General Information

1.3.1. Scope of Delivery

The delivery includes the following items:



a. Spin Air Basic main unit (Cat. 10005532)
b. Power supply (Cat. 90001566) and power cord.
c. Allen wrench (Cat. 90005516)

Apart from the list of items above, the delivery also comes with the Spin Air Basic Datamatrix (Doc. 50008875).

See APPENDIX A Ordering Codes for more information.
NOTE: This is an image of the standard model. The set head of the Spin Air Basic depends on the order of the customer and the appearance.

1.3.2. Technical Assistance

At IUL we pride ourselves on the quality and availability of our technical support. Our Technical Service Department is staffed by experienced technicians with extensive practical and theoretical expertise in the use of IUL products. If you have any questions or experience any difficulties regarding the Spin Air Basic or IUL products in general, do not hesitate to contact us or your local distributors.

IUL customers are a major source of information regarding advanced or specialized uses of our products. This information is helpful to other customers as well as to the researchers at IUL. We, therefore, encourage you to contact us if you have any suggestions about product performance or new applications and techniques.

For technical assistance, contact the IUL Technical Service Department or your local distributors.

1.3.3. Policy Statement

It is the policy of IUL to improve products as new techniques and components become available. IUL reserves the right to change the specifications of products at any time.

To produce useful and appropriate documentation, we appreciate your comments on these Instructions for Use. Please contact IUL Technical Service with any feedback.

1.4. Requirements for Spin Air Basic Users

Table 1 covers the general level of competence for the use and servicing of the Spin Air Basic.

Task	Personnel	Training and experience
Routine use	Laboratory technicians or equivalent	Trained in techniques for laboratory instrument operation
Servicing	IUL Service Specialists only	Trained, certified, and authorized by IUL

2. SAFETY INFORMATION

Before using the Spin Air Basic, it is essential that you read these Instructions for Use carefully. Following the instructions and safety information in this Instructions for Use will ensure safe operation and maintain the system in a safe condition.

The following types of safety information appear throughout the Spin Air Basic Instructions for Use.

WARNING 	The term WARNING is used to inform you about situations that could result in **personal injury** to you or other persons. Details about these circumstances are given in a box like this one.

CAUTION 	The term CAUTION is used to inform you about situations that could result in **damage to the instrument** or other equipment. Details about these circumstances are given in a box like this one.

The advice given in this Instructions for Use is intended to supplement, not supersede, the normal safety requirements prevailing in the user's country.

2.1. Proper Use

WARNING/ CAUTION 	**Risk of personal injury and material damage** Improper use of the Spin Air Basic instrument may cause personal injury or damage to the instrument. The instrument must only be operated by qualified personnel.

CAUTION 	**Damage to the instrument** Avoid spilling water or chemicals onto the Spin Air Basic instrument. Damage caused by water or chemical spillage will void your warranty.

In case of emergency, switch off the Spin Air Basic and unplug the power cord from the power outlet.

2.2. Electrical Safety

If the operation of the Spin Air Basic is interrupted in any way (e.g., due to interruption of the power supply or a mechanical error), first switch off the instrument using the power button, then unplug the power cord from the power outlet. Contact IUL or IUL's Partner Technical Service after such an incident.

WARNING 	**Rechargeable Batteries** This WARNING is applicable to routine users as well as to technical service users. • Do not dismantle, open or shred batteries. • Keep batteries out of the reach of children. • Seek medical advice immediately if a battery has been swallowed. • Do not expose batteries to heat or fire. Avoid storage in direct sunlight. • Do not remove a battery from its original packaging. • Do not subject batteries to mechanical shock • In event of a cell leaking, do not allow the liquid to encounter the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice. • Only use Batteries provided by IUL or your local distributor. After extended periods of storage, it may be necessary to charge and discharge the batteries several times to obtain maximum performance.

WARNING 	**Electrical hazard** Any interruption of the protective conductor (earth/ground lead) inside or outside the instrument or disconnection of the protective conductor terminal is likely to make the instrument unsafe. **This must be checked after service or maintenance.** Intentional interruption is prohibited. **Lethal voltages inside the instrument** This WARNING applies to routine users as well as to technical service users. • Risk of electrical shock and energy hazard. All failures should be examined by a qualified technician. Please do not remove the case of the AC adaptor by yourself! • Adaptors should be placed on a reliable surface. A drop or fall could cause damage. • Please do not place the AC adaptor in places with high moisture or near the water. • Please do not place the AC adaptor in places with high ambient temperature or near the fire source. About the maximum ambient temperature, please refer to "APPENDIX B Technical Data".

	• Disconnect the AC adaptor from the AC power before cleaning. Do not use any liquid or an aerosol cleaner. Only use a damp cloth to wipe it. In case of replacement or a loosing of the AC adaptor or the mains power cord, these must be replaced only with the AC adaptors or power cords listed on the "Ordering Codes" In case of replacement, this must be ordered to IUL's Partner.

To ensure satisfactory and safe operation of the Spin Air Basic:

- The power cord of the external power supply must be connected to a line power outlet that has a protective conductor (earth/ground).
- No other external power supply neither power cords than the specified in "APPENDIX A Ordering Codes" must be used. In case of replacement, this must be ordered to IUL.
- The instrument must not be operated with the head removed.
- If you suspect any instrument damage, contact IUL Technical Services or your local distributor.

If the Spin Air Basic becomes electrically unsafe, prevent other personnel from operating it, and contact IUL Technical Service.

The instrument may be electrically unsafe if:

- The instrument or the external power supply is shown to be damaged.
- The instrument has been stored under unfavourable conditions for a prolonged period.
- A different external power supply is used other than the one provided by IUL.

WARNING 	**Risk of electric shock** In case of replacement or loss of the external power supply or the power cord, these must be replaced only with the external power supply or power cords listed on the "APPENDIX A Ordering Codes" and provided by IUL.

2.3. Environment

Operating conditions

CAUTION 	**Explosive atmosphere** The Spin Air Basic is not designed for use in an explosive atmosphere.
CAUTION 	**Overheating risk** Do not cover gaps and slots that ensure the ventilation of the Spin Air Basic.

2.4. Biological Safety

Use safe laboratory procedures as outlined in publications such as Biosafety in Microbiological and Biomedical Laboratories, HHS:

https://www.cdc.gov/labs/pdf/CDC-BiosafetyMicrobiologicalBiomedicalLaboratories-2009-P.PDF

See: Section III—Principles of Biosafety

Laboratory Practices and Technique

The most important element of containment is strict adherence to standard microbiological practices and techniques. Persons working with infectious agents or potentially infected materials must be aware of potential hazards and must be trained and proficient in the practices and techniques required for handling such material safely. The director or person in charge of the laboratory is responsible for providing or arranging the appropriate training of personnel.

WARNING 	**Samples containing infectious agents** Some samples used with the Spin Air Basic may contain infectious agents. Handle such samples following the required safety regulations. The responsible person(s) (e.g., laboratory manager) must take the necessary precautions to ensure that the workplace is safe and that the instrument operators are suitably trained and not exposed to hazardous levels of infectious agents, as defined in the applicable Safety Data Sheets (SDSs) or OSHA[1] ACGIH[2] or COSHH[3] documents. Venting of fumes and disposal of wastes must be following all national, state, and local health and safety regulations and laws.

2.5. Chemicals

WARNING 	**Hazardous chemicals** Some chemicals used with the Spin Air Basic may be hazardous. Always wear safety glasses, gloves, and a lab coat. The responsible person(s) (e.g., laboratory manager) must take the necessary precautions to ensure that the workplace is safe and that the instrument operators are suitably trained and not exposed to hazardous levels of toxic substances (chemical or biological), as defined in the applicable Safety Data Sheets (SDSs) or OSHA, ACGIH or COSHH documents. Venting of fumes and disposal of wastes must be following all national, state, and local health and safety regulations and laws.

2.6. Maintenance Safety

CAUTION 	**Damage to the instrument** Do not use spray bottles containing alcohol or disinfectant to clean the surface of the Spin Air Basic instrument. Do not use products containing alcohol or other corrosive solvents to clean the Spin Air Basic instrument.
WARNING 	**Risk of electric shock** Do not open the panels on the instruments. **Risk of personal injury and material damage** Only perform maintenance that is specifically described in this Instructions for Use.

[1] OSHA: Occupational Safety and Health Administration (United States of America).
[2] ACGIH: American Conference of Government Industrial Hygienists (United States of America).
[3] COSHH: Control of Substances Hazardous to Health (United Kingdom).

2.7. Waste Disposal

Used consumables, such as sample dishes, may contain hazardous chemicals or infectious agents. Such waste must be collected and disposed of properly following local safety regulations.

For disposal of waste electrical and electronic equipment (WEEE) see “APPENDIX D Waste Electrical and Electronic Equipment (WEEE)”.

2.8. Symbols on the Spin Air Basic Device

Symbol	Location	Description
	Type plate on the bottom of the device.	Waste Electrical and Electronic Equipment (WEEE), see “APPENDIX D Waste Electrical and Electronic Equipment (WEEE)”.
CE	Type plate on the bottom of the device.	APPENDIX I EU Declaration of Conformity.

3. GENERAL DESCRIPTION

3.1. Device Overview

Front view:

A. Head


B. Display

C. Keyboard



Rear view:

D. Holder and handle

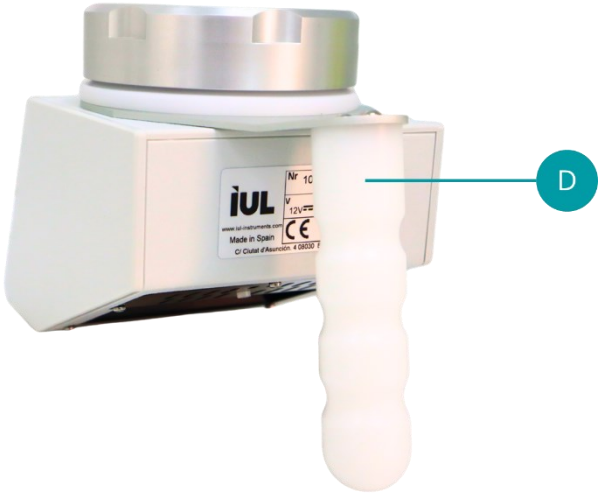


Bottom view:

E. Air exhaust

F. Power supply connector

G. Reset buttonhole

H. Coupling tripod with base

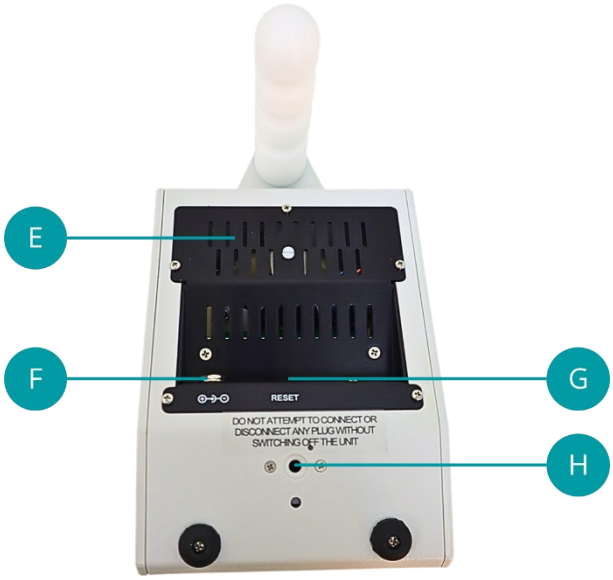

See “APPENDIX A Ordering Codes” for more information.

NOTE: Use only the AC Adapter supplied by IUL to charge the Spin Air Basic.

When the device shows a “blank screen” and it is not possible to switch on the device either start the charging process, do the following: Connect the charger, wait 10 seconds, and then push the “Reset” button. If the problem isn’t solved, contact your Technical Service.

3.2. Control Panel

KEY	Action		
	Pressed **once**	Pressed longer than **3 sec**	Pressed **continuously**
	START a program	Switch ON the instrument if it is off	NA
	STOP running a program RETURN TO PREVIOUS menu	OFF	NA
	Move ONE STEP forward READ BARCODE if the Barcode Scanner is connected	NA	STEP FORWARD quickly
	CHANGE menu ACCEPT parameter	NA	NA
	ONE STEP BACK	NA	STEP BACK quickly

A sample process can be always interrupted by pushing  . The icon showing the time on the left will start flashing. The interrupted process cannot be restarted. Pushing any key will return to the main menu.

3.3. Accessories of the Device

3.3.1. AC Adaptor

The AC Adaptor (Cat. 90001566) is included. It converts the Mains AC supply to a DC low voltage supply required by the Spin Air Basic.

3.3.2. Power Cord

The Spin Air Basic is equipped with a power cord with a plug suitable for the destination country. See “APPENDIX A Ordering Codes” for details.

CAUTION 	**Risk of malfunction or harmful interference** Conformity to **2014/30/EU** and FCC rules could be compromised by the use of an AC adaptor or power cord not provided by IUL see “APPENDIX A Ordering Codes” and “APPENDIX E RoHS Statement”.

WARNING 	**Risk of electric shock** In case of replacement or a loosing of the power cord, this must be replaced only with one of the power cords listed on the “APPENDIX A Ordering Codes” delivered by IUL.

4. INSTALLATION

4.1. Unpack the Device

The packaging of the Spin Air Basic can be stored for reuse.





4.2. Site Requirements

Place the Spin Air Basic device on a stable surface, far away from powerful lights and near an earthed/grounded electrical outlet. If a surface is not available, connect the Spin Air Basic to a tripod. See “APPENDIX A Ordering Codes”.

For a long period of use, kept the Spin Air Basic device connected to the power supply.

4.3. Power Cable Connection

The socket for connecting the power cable is on the back of the Spin Air Basic device.

When the Spin Air Basic is not in use for a long period, we recommend disconnecting the power cable.


5. OPERATING PROCEDURES

5.1. Sampling

1. witch on the instrument by pressing the button for more than 3 seconds. The software version will be shown:



In stand-by mode, the following display indications will be shown:

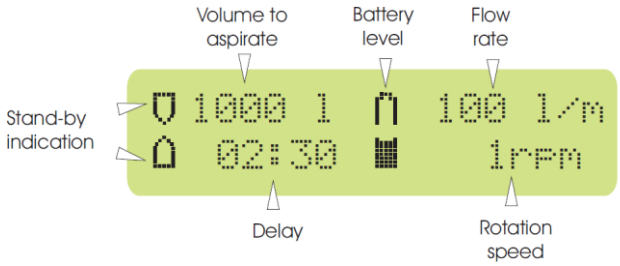

2. Select the parameters before running the program by pressing the arrow up ( ), down arrow ( ), and tick keys ( ):


3. Adjust the volume (10 to 9900L).



4. Adjust the start delay in minutes and seconds.




5. Adjust the rotation speed (0 to 4 rpm).

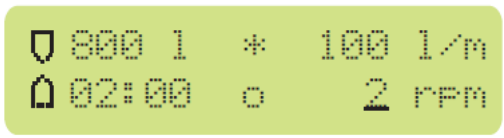


6. Ensure that the current air flow corresponds to the head used. Use 100l/m for 90mm set head and 60l/min for 60mm set head.

7. Run the program by pressing :

At the start of the operation, all the “sand” is displayed in the upper part of the hourglass on the left-hand screen. As sampling progresses, the “sand” moves to the lower part in proportion to the volume of air sampled.

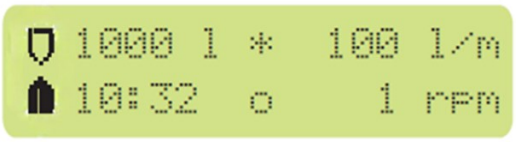

The remaining and the elapsed time will be indicated as is shown below:

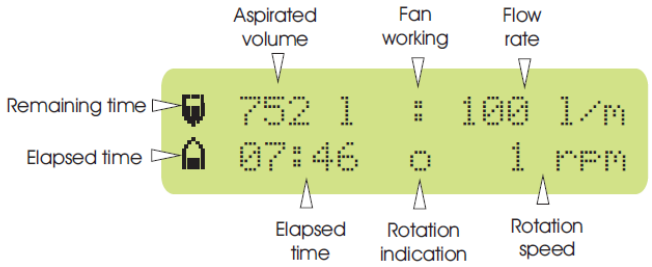


8. Switch off the device by pressing the cross key for more than 3 seconds. The device will be automatically switched off if it is not used for more than 5 minutes.


NOTE: A sample process can always be interrupted by pressing the button. The hourglass icon will start flashing with last state. The interrupted process cannot be restarted. By pressing any key, it will return to the main menu.

5.2. Battery Charging


The instrument can be operated while connected to the charger or while running on battery power. During charging, the sampling parameters remain visible, and the battery indicator displays the charging status. The battery is fully charged after approximately 3 hours.

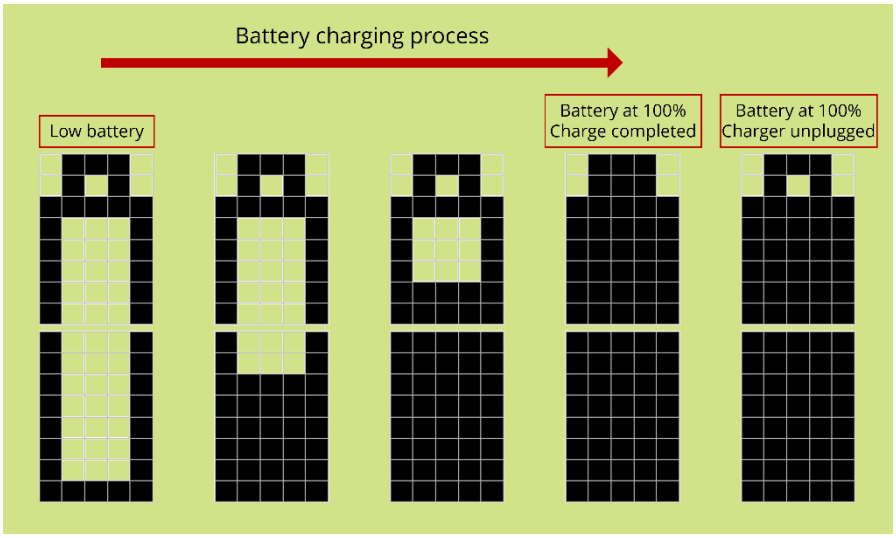

When the battery is too low, the device will display the message “Battery Low” and the battery icon will blink, as shown in the following figure:

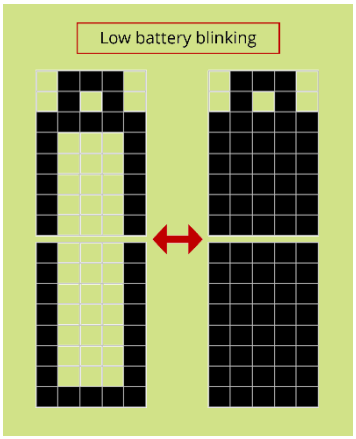


If a problem occurs during the battery charging process, the battery icon will blink, as shown in the figure below. In this case, check the power supply. If the power supply is functioning correctly and the battery icon continues to blink, contact Technical Support.

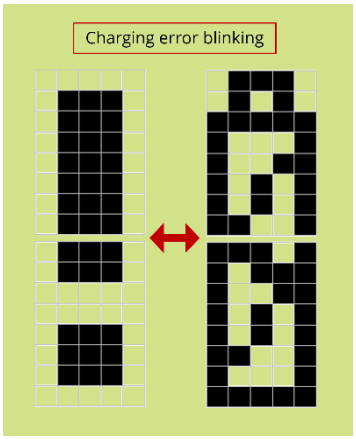

5.3. Advanced Menu

Push the tick key ( ) for more than 2 seconds to access the advanced menu:

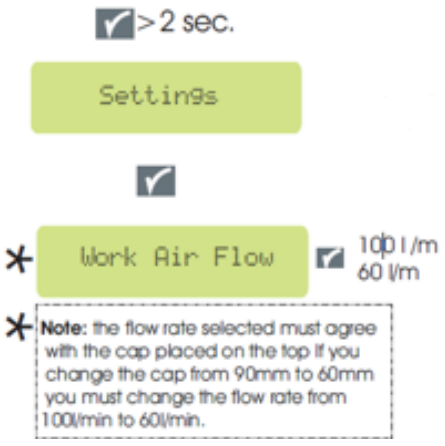

To return one step back push the cross key ( ).

6. CALIBRATION

The instrument calibration is valid for 2 years or 1 million air liters. When the instrument reaches one of these situations it shows a message to every operation notifying about that.

The calibration must be performed by a trained operator, and it is compulsory to use an anemometer manufactured for this specific purpose. Contact your IUL distributor when calibration is needed.

7. TROUBLESHOOTING

Symptom	Probable cause	Recommended action
The unit cannot turn ON.	Battery discharged.	Plug the power supply unit. Wait 5 seconds. If the battery icon doesn't appear then press the RESET button.
The unit switches off spontaneously.	The mains cord is unplugged and the battery getting too low.	Plug the mains supply and let the battery fully charge.
"Power flow Err" is displayed.	The Dish is obstructed / The battery voltage has dropped / No head.	Clean the Dish holes, dry them and place properly the Dish on the device / Check the charging process and if you believe a part is damaged or the battery is in bad condition contact your distributor or the Technical Service / Verify the installation of the head.
"Calibration Needed" is displayed.	The calibration is out of limits (2 years or 1 million liters).	Perform the calibration process.
"Battery Low" is displayed.	The battery is empty.	Plug the charger and wait about 4 hours.

8. MAINTENANCE

8.1. Cleaning Procedure

To clean the Spin Air Basic device use a damp brush or cloth moistened with a solution of hot water and soap. The stainless steel and aluminum heads can be autoclavable or clean with an alcohol solution (70%).

CAUTION 	**Damage to the device.** Do not sprinkle inside the device.

APPENDIX A ORDERING CODES

Spare Parts

Product	Contents	Cat. #
Power cord (UK)	1 unit	90002497
Power cord (Argentina)	1 unit	90002499
Power cord (EU/Schuko)	1 unit	90002495
Power cord (China/Australia)	1 unit	90002498
Power cord (Japan/USA)	1 unit	90002496
Adapter AC	1 unit	90001566

Optional Accessories

Product	Contents	Cat. #
Aluminium head and lid for 90 mm Petri dishes	1 unit	90005504
Aluminium head and lid for 60 mm Petri dishes (RODAC)	1 unit	90005505
Calibration kit with transport case	1 unit	900012093
Tripod with base	1 unit	90005511
Plastic head and lid for 90 mm Petri dishes	1 unit	90005525
Plastic head and lid for 60 mm Petri dishes (RODAC)	1 unit	90005526
INOX head and lid for 90 mm Petri dishes	1 unit	90005527
Tripod base	1 unit	900011099

APPENDIX B TECHNICAL DATA

IUL reserves the right to change specifications at any time.

Specification

Airflow	100 L/min (90 mm plate) 60 L/min (55–60 mm and RODAC plate) Controlled by microprocessor
Air Total Volume	10 – 9900 L
Total Air Volume Reference	The total air volume parameter must be referenced to sea-level pressure
Delay to Start	Up to 60 minutes, adjustable in seconds
Rotation Speed	0, 1, 2, 3, and 4 rpm
Tripod Thread	Located at the bottom
Communication	RS-232C to printer / DB9M–USB to PC
Data Export format	XML, CSV
Barcode Connection	Compatible with IUL Barcode Reader
IUL Switching Adapter (90001566)	Input: 100–240 V AC, 50/60 Hz ; Output: 12 V DC, 15 W
Battery Pack	Nickel-metal hydride, 7.2 V
Battery charging time	Approx. 3h
Range	8 hours full charge
Automatic power-off	After 5 min of inactivity

Operating Conditions

Power Range	12 V ±3% Supplied with an IUL external power supply unit suitable for 100–240 V
Frequency Range	DC Supplied with an external power supply unit suitable for 50–60 Hz
Maximum Power	18 W
Overvoltage Category	II
Air Temperature	10–40 °C / 50–104 °F

Relative Humidity	10–75% (non-condensing)
Atmospheric Pressure	1 atm ±5% (101.3 kPa)
Place of Operation	Indoor use only
Pollution Level	2
Energy Efficiency Level	VI

Transportation conditions

Air Temperature	10–40 °C / 50–104 °F, in the manufacturer's packaging
Relative Humidity	Maximum 75% (non-condensing)

Storage conditions

Air Temperature	10–40 °C / 50–104 °F, in the manufacturer's packaging
Relative Humidity	Maximum 75% (non-condensing)

Weight and dimensions

Dimensions (W×H×D)	145 × 190 × 215 mm / 5.71 × 7.48 × 8.46 in
Weight	Spin Air basic: 1.7 kg / 3.75 lb Carrying case: 2.24 kg / 4.85 lb (empty)

APPENDIX C WARRANTY

The IUL instrument has a 12-month warranty period.

The warranty is void when misuse of the equipment can be proved. Damage or faults caused by impacts, chemical or corrosive products, liquids, dampness, or other external factors, such as radiation, fire, or inadequate transport are not included.

In addition, the warranty will not apply if the equipment has been handled, repaired, or modified by unqualified or specifically designated personnel.

As part of this warranty, all shipping costs to the selected technical service are the responsibility of the Buyer.

APPENDIX D WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE)

This section provides information about the disposal of waste electrical and electronic equipment by users.

The crossed-out wheeled bin symbol (see below) indicates that this product must not be disposed of with other waste. The product must be disposed to a certified treatment facility or a recycling collection point, according to local legislation.

The separate collection and recycling of electronic waste equipment at the time of disposal helps conserve natural resources and ensures that the product is recycled, protecting human health and the environment.



APPENDIX E ROHS STATEMENT

The following information has been made available to comply with The Restriction of Hazardous Substances Directive, (RoHS2 & RoHS3), short for Directive on the restriction of the use of certain hazardous substances in electrical and electronic equipment.



APPENDIX F STATISTICAL CORRECTIONS WHEN SAMPLING WITHOUT DISK ROTATION

The principle of Spin Technology is mainly based on two features: the rotating Petri dish and the distribution of the holes. The rotating Petri dish feature provides a genuine count without the need to apply statistical corrections. This is because a motionless system only uses 5% of the total surface of the agar. The IUL's Spin Air has a unique rotation motion device combined with a distinctive distribution of the holes that covers 100% of the agar surface.

If the Spin Air is used **without rotation** (0 rpm) a correction of the result must be applied. This correction takes into account the probability that one microorganism travels across the same hole that a previous one. This correction that is not important for low counts becomes dramatically high when there are more than 150 CFU per plate.

The statistical treatment due to Feller (1950) says that the probable statistically total (Pr) can be calculated according to the number of holed (N) and the real count (r).

$Pr = N \cdot [1/N + 1/(N-1) + 1/(N-2) + .. 1/(N – r + 1)]$

Applying the previous formula to the Spin Air for 90 mm plate the probable total count (Tc) is when (C) is the real count.

$Tc = 400 \cdot [1/400 + 1/399 + 1/398 + ... + 1/(400 – C + 1)]$

Up to 20 colonies, there is no difference between the motionless and the rotation system. At 50 counts you must add around 5%. At 70 the difference is 10%. If the count reaches 100 CFUs you must add 15% of the total count. This percentage reaches 25% at 150 counts and it's more than 100% at more than 300 counts.


Using the **rotation motion** of Spin Air all these calculations must be avoided.

APPENDIX G TABLE OF STATISTICAL CORRECTIONS FOR SAMPLING USING A 55/60 MM DISH

(USE ONLY IN CASE SAMPLES TAKEN WITHOUT ROTATION)

Total number of holes: N = 240
Colony number counter: r
Colony number estimated: PR

r	Pr	r	Pr	r	Pr	r	Pr	r	Pr
1	1	51	57	101	131	151	237	201	434
2	2	52	58	102	132	152	240	202	440
3	3	53	60	103	134	153	243	203	446
4	4	54	61	104	136	154	245	204	452
5	5	55	62	105	138	155	248	205	459
6	6	56	64	106	139	156	251	206	466
7	7	57	65	107	141	157	254	207	473
8	8	58	66	108	143	158	257	208	480
9	9	59	68	109	145	159	260	209	488
10	10	60	69	110	147	160	263	210	496
11	11	61	70	111	149	161	266	211	504
12	12	62	72	112	150	162	269	212	512
13	13	63	73	113	152	163	272	213	520
14	14	64	74	114	154	164	275	214	529
15	15	65	76	115	156	165	278	215	539
16	17	66	77	116	158	166	281	216	548
17	18	67	78	117	160	167	285	217	558
18	19	68	80	118	162	168	288	218	569
19	20	69	81	119	164	169	291	219	579
20	21	70	83	120	166	170	295	220	591
21	22	71	84	121	168	171	298	221	603
22	23	72	85	122	170	172	301	222	616
23	24	73	87	123	172	173	305	223	629
24	25	74	88	124	174	174	309	224	643
25	26	75	90	125	176	175	312	225	658
26	27	76	91	126	178	176	316	226	674
27	29	77	93	127	180	177	320	227	691
28	30	78	94	128	182	178	323	228	710
29	31	79	96	129	184	179	327	229	730
30	32	80	97	130	187	180	331	230	751
31	33	81	99	131	189	181	335	231	775
32	34	82	100	132	191	182	339	232	802
33	35	83	102	133	193	183	343	233	832
34	37	84	103	134	195	184	348	234	866
35	38	85	105	135	198	185	352	235	906
36	39	86	106	136	200	186	356	236	954
37	40	87	108	137	202	187	361	237	1014
38	41	88	109	138	205	188	365	238	1094
39	42	89	111	139	207	189	370	239	1214
40	44	90	113	140	209	190	375	240	1454
41	45	91	114	141	212	191	379		
42	46	92	116	142	214	192	384		
43	47	93	117	143	217	193	389		
44	48	94	119	144	219	194	394		
45	50	95	121	145	222	195	400		
46	51	96	122	146	224	196	405		
47	52	97	124	147	227	197	410		
48	53	98	126	148	229	198	416		
49	55	99	127	149	232	199	422		
50	56	100	129	150	235	200	428		

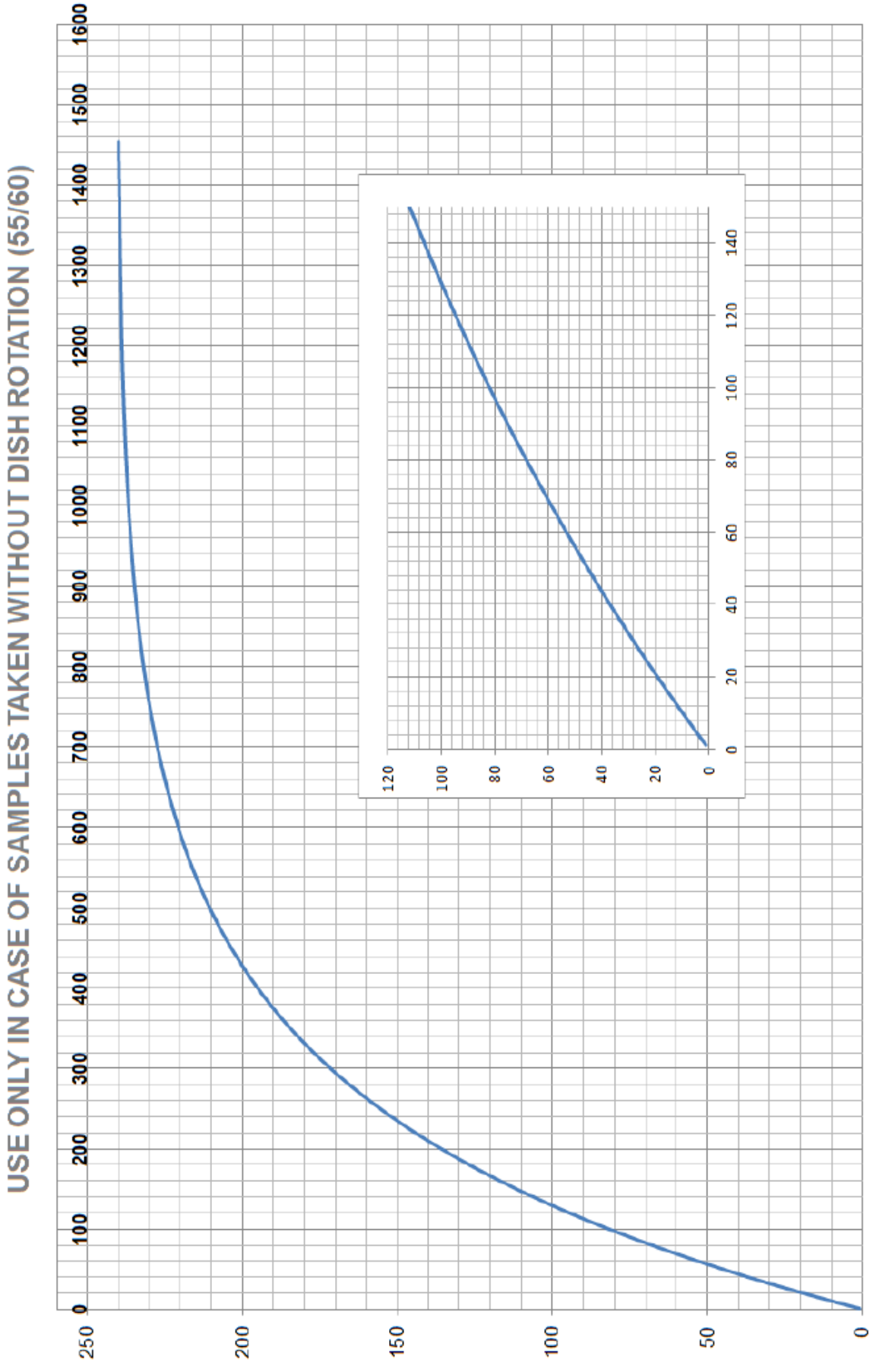
USE ONLY IN CASE OF SAMPLES TAKEN WITHOUT DISH ROTATION (55/60)
0
100
200
300
400
500
600
700
800
900
1000
1100
1200
1300
1400
1500
1600
250
200
150
100
50
0
0
20
40
60
80
100
120
140
120
100
80
60
40
20
0

APPENDIX H TABLE OF STATISTICAL CORRECTIONS FOR SAMPLING USING A 90 MM DISH


(USE ONLY IN CASE SAMPLES TAKEN WITHOUT DISH ROTATION)

Total number of holes: N = 400

Colony number counter: r

Colony number estimated: PR

r	Pr	r	Pr	r	Pr	r	Pr	r	Pr	r	Pr	r	Pr	r	Pr
1	1	51	54	101	116	151	189	201	279	251	394	301	557	351	836
2	2	52	56	102	118	152	191	202	281	252	397	302	561	352	844
3	3	53	57	103	119	153	193	203	283	253	400	303	565	353	853
4	4	54	58	104	120	154	194	204	285	254	402	304	569	354	861
5	5	55	59	105	122	155	196	205	287	255	405	305	573	355	870
6	6	56	60	106	123	156	197	206	289	256	408	306	578	356	879
7	7	57	61	107	124	157	199	207	291	257	411	307	582	357	888
8	8	58	63	108	126	158	201	208	293	258	413	308	586	358	897
9	9	59	64	109	127	159	202	209	295	259	416	309	591	359	907
10	10	60	65	110	128	160	204	210	297	260	419	310	595	360	917
11	11	61	66	111	130	161	206	211	299	261	422	311	599	361	927
12	12	62	67	112	131	162	207	212	301	262	425	312	604	362	937
13	13	63	68	113	133	163	209	213	304	263	428	313	608	363	947
14	14	64	70	114	134	164	211	214	306	264	431	314	613	364	958
15	15	65	71	115	135	165	212	215	308	265	433	315	618	365	969
16	16	66	72	116	137	166	214	216	310	266	436	316	622	366	981
17	17	67	73	117	138	167	216	217	312	267	439	317	627	367	992
18	18	68	74	118	140	168	218	218	314	268	442	318	632	368	1005
19	19	69	76	119	141	169	219	219	317	269	445	319	637	369	1017
20	20	70	77	120	142	170	221	220	319	270	449	320	642	370	1030
21	22	71	78	121	144	171	223	221	321	271	452	321	647	371	1043
22	23	72	79	122	145	172	224	222	323	272	455	322	652	372	1057
23	24	73	80	123	147	173	226	223	325	273	458	323	657	373	1071
24	25	74	82	124	148	174	228	224	328	274	461	324	662	374	1086
25	26	75	83	125	150	175	230	225	330	275	464	325	667	375	1102
26	27	76	84	126	151	176	232	226	332	276	467	326	673	376	1118
27	28	77	85	127	153	177	233	227	335	277	471	327	678	377	1134
28	29	78	87	128	154	178	235	228	337	278	474	328	684	378	1152
29	30	79	88	129	156	179	237	229	339	279	477	329	689	379	1170
30	31	80	89	130	157	180	239	230	342	280	480	330	695	380	1189
31	32	81	90	131	158	181	241	231	344	281	484	331	701	381	1209
32	33	82	92	132	160	182	242	232	346	282	487	332	706	382	1230
33	34	83	93	133	161	183	244	233	349	283	491	333	712	383	1252
34	35	84	94	134	163	184	246	234	351	284	494	334	718	384	1276
35	37	85	95	135	164	185	248	235	353	285	497	335	724	385	1301
36	38	86	97	136	166	186	250	236	356	286	501	336	730	386	1327
37	39	87	98	137	167	187	252	237	358	287	504	337	737	387	1356
38	40	88	99	138	169	188	254	238	361	288	508	338	743	388	1387
39	41	89	101	139	171	189	255	239	363	289	511	339	749	389	1420
40	42	90	102	140	172	190	257	240	366	290	515	340	756	390	1456
41	43	91	103	141	174	191	259	241	368	291	519	341	763	391	1496
42	44	92	104	142	175	192	261	242	371	292	522	342	769	392	1541
43	45	93	106	143	177	193	263	243	373	293	526	343	776	393	1591
44	47	94	107	144	178	194	265	244	376	294	530	344	783	394	1648
45	48	95	108	145	180	195	267	245	378	295	534	345	791	395	1715
46	49	96	110	146	181	196	269	246	381	296	537	346	798	396	1795
47	50	97	111	147	183	197	271	247	384	297	541	347	805	397	1895
48	51	98	112	148	185	198	273	248	386	298	545	348	813	398	2028
49	52	99	114	149	186	199	275	249	389	299	549	349	820	399	2228
50	53	100	115	150	188	200	277	250	391	300	553	350	828	400	2628

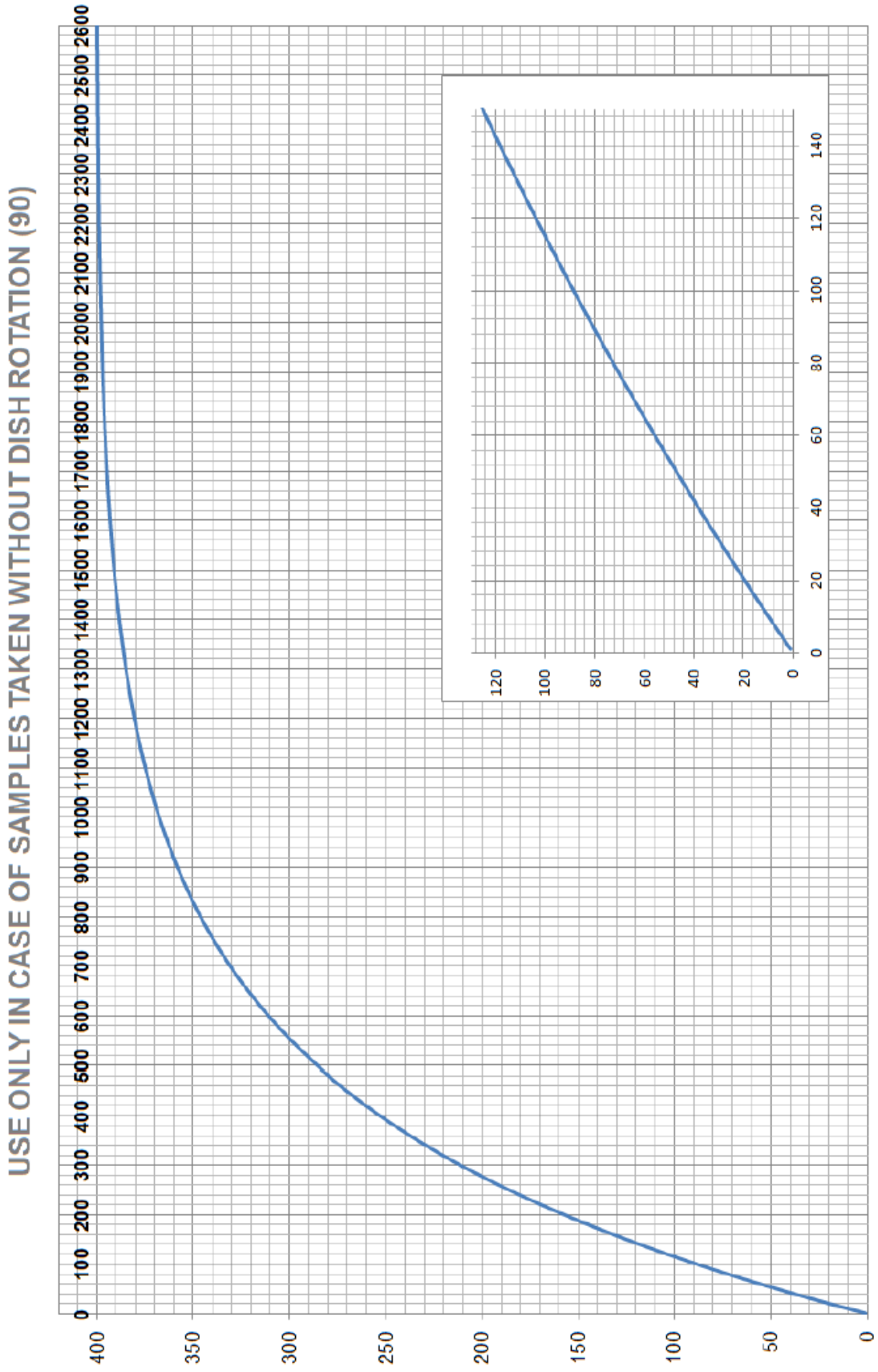
USE ONLY IN CASE OF SAMPLES TAKEN WITHOUT DISH ROTATION (90)
0 100 200 300 400 600 600 700 800 900 1000 1100 1200 1300 1400 1500 1600 1700 1800 1900 2000 2100 2200 2300 2400 2500 2600
400 350 300 250 200 150 100 50 0
0 20 40 60 80 100 120 140
120 100 80 60 40 20 0

APPENDIX I EU DECLARATION OF CONFORMITY

Manufacturer Name:
IUL
Manufacturer Address:
Carrer de la Ciutat d'Asuncion, 4 08030 Barcelona Spain

Declares that the product Name:
Spin Air Basic

Catalog Number:
90005500
Conforms the Council Directives:

- Directive **2014/30/EU** of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility (recast).
- Directive **2014/35/EU** of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits (recast).
- Directive **2011/65/EU** of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment and its amendments.

Standards which Conformity is Declared:

- **EN 61326-1:2013** Electrical equipment for measurement, control and laboratory use - EMC requirements -- Part 1: General requirements (IEC 61326-1:2012)
- **EN 61010-1:2010** Safety requirements for electrical equipment for measurement, control, and laboratory use -- Part 1: General requirements.

Barcelona, January the 23rd 2024

CE

Daniel Blanch / Director of Electronical Engineering

APPENDIX J CERTIFICATE OF COMPLIANCE EN 17141:2020

UK Health Security Agency

Certificate of Compliance Testing

Biological Efficiency Testing of the IUL Spin Air 100 and 60 Samplers Using Techniques Described in EN 17141:2020 Annex E, Section E6.2

Report No. 23/011

Report Prepared For: IUL Instruments

Issue Date: 10 January 2024

Test Summary

Tests were undertaken using the simplified laboratory method (Annex E6.2) sampling environmental aerosols. Compared to the reference slit sampler the IUL Spin Air 100 was found to be 83.30% efficient and the Spin Air 60 was found to be 86.41% efficient. This shows compliance to acceptance criteria detailed in Annex E.6.4. Tests were also undertaken using a variation of the aerosol chamber method, comparing the ratio of the recoveries of mixed aerosol of *Staphylococcus epidermidis* and *Bacillus atrophaeus* spores. It was found that the compared to the reference sampler the Spin Air 100 was 91.73% efficient and the Spin Air 60 was 98.49% efficient.

Report Written By

Name: Mr Simon Parks

Title: Senior Biosafety Scientist

Report Authorised By

Name: Ms Helen Hookway

Title: Senior Biosafety Scientist

Please be aware that the use of the Royal Coat of Arms is highly restricted and cannot be copied. Please do not put the UKHSA logo on your website or use our name to endorse your products. Any reference to UKHSA needs to be approved by us before it can be used.

IUL TECHNICAL SUPPORT CONTACT

Global HQ	
IUL	
Phone:	+34 932 740 232 ext: 4318
Email:	support@iul-inst.com
Address:	C/ Ciutat d'Asunción, 4 08030 Barcelona Spain

United States Operations	
IUL Inc.	
Phone:	+1 (571) 237-0188
Email:	service@iul-instruments.us
Address:	4700 Eisenhower Ave, ALX-F2 Alexandria, VA, 22304 USA

Germany, Austria, the Netherlands Operations	
IUL GmbH	
Phone:	+49 7152 615 0998
Email:	service.de@iul-inst.com
Address:	Hertichstr. 10 71229 Leonberg Germany